


Powered by



Hemo Vet

Hemoglobin testing for veterinary care

Scan here to view

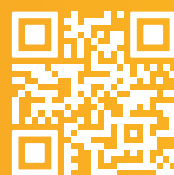

product webpage

Hemoglobin and hematocrit results in one test

Evaluated for various species

Reliable in a wide temperature range

Adjustable display for different light conditions

EKF
Diagnostics

Hemo Vet

Easy-to-use

- User-friendly features minimize training time
- Step-by-step instructions on screen
- Adjustable backlit touch screen - readable in intense sunlight or bright artificial light
- User-selectable language menu
- Soft-load cuvette holder eliminates risk of contamination

Practical and portable

- Hemoglobin and hematocrit results from one sample available in 25-60 seconds (depending on the concentration)
- Venous, arterial or capillary blood
- Sample volume only 8µl depending on cuvette type
- Microcuvette minimises air bubbles
- Compact and light (700g)
- No maintenance required. Auto self-test
- Integrated rechargeable battery (100 hours)

Accurate and reliable

- Operating ambient temperature 15°C - 40°C
- Photometric, azide methemoglobin method
- Measuring range: 0 - 25.6g/dL; 0 - 15.9mmol/L
- Imprecision: CV<2%
- Linearity: 0-20g/dl ± 0.3g/dl; >20g/dl ± 0.7g/dl
- Calibrated to NCCLS reference method
- Stores up to 4,000 results
- Control cuvette included

Catalog no.


3032-0010-0219

0760-0102

0760-0104

Name

- Hemo Vet
- EKF Link license
- EKF Link Lite

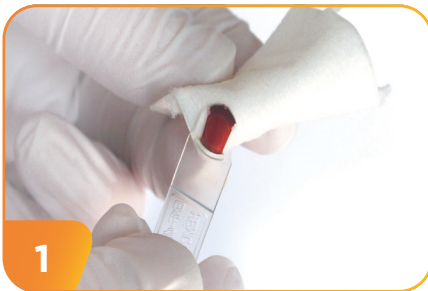


1

Collect blood sample.



2

Put microcuvette into analyser.

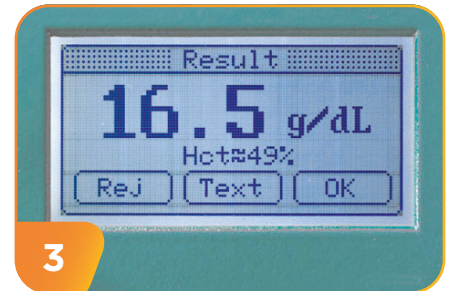

3

Result appears in 25 - 60 seconds.



Manufacturer

EKF-diagnostics GmbH
Ebendorfer Chaussee 3
39179 Barleben Germany
Tel: +49 (0) 39203 511 0
Email: sales@ekfdiagnostics.com
ekfdiagnostics.com

Distributed by